


FagronLab™ LF-1300 UV

Positive pressure workbench

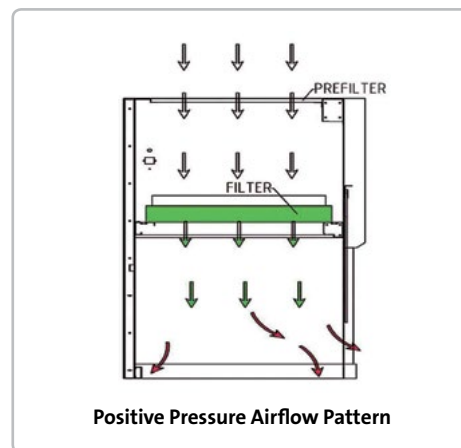

FagronLab™ LF-1300 UV

- New Air Purifier System necessary for compounding and homeopathic laboratories
- Through a special HEPA filtration system with efficiency of 99.995 % at 0,3 um, it flows unidirectional, filtered streams of air with a steady uniform flow to the working area
- Maintains positive air pressure protecting your preparations from invading particles and the intrusion of contaminated room air
- Provides an excellent environment for compounding and homeopathic forms with high safety demands
- Antibacterial protection with UV lamp. Easy to use operational digital LCD display
- Working surface made of high quality stainless steel, without gaps in joints to prevent bacteria accumulation. Internal lamp for lighting compounding area
- Automatically motorized front window
- Wooden table with 6 drawers and special position for a sterilization furnace (Sterilization furnace and equipment shown on the picture - extra)
- Warranty 1 year

FagronLab™ LF-1300 UV Homeo edition

External Size (HxLxD): 236x137,5x75cm

Internal Size (HxLxD): 61x120x65cm

Weight: 294kg

FagronLab™ LF-1300 UV Pharma edition

External Size (HxLxD): 236x180x70cm

Internal Size (HxLxD): 61x170x65cm

Weight: 379kg